



Pacific imports manufacturing inc.PIM

## **Disposable medical (three-layer) mask**

### **Product markings and instructions**

Pacific Import Manufacturing, Inc.  
13000 Temple Avenue,  
City of Industry, CA 91746

# Attention!

Disposable medical (three-layer) mask

Disposable protective (three-layer) mask

No difference in materia

**The biggest difference is that before making a disposable medical (three-layer) mask, it needs to be controlled and must be disinfected in the disinfection chamber for 14 days before it can be used to make a qualified medical mask**

# Disposable medical (three-layer) mask

**This product implements GB / T32610-2016 standard, and its main performance is as**

**follows:**

- (1) The filtration efficiency of non-oily particles is  $\geq 90\%$ , and the filtration efficiency of oily particles is  $\geq 80\%$ .
- (2) Protective effect  $\geq 65\%$ .
- (3) The breaking strength at the connection point between the earband and the protective body is  $\geq 20\text{N}$ .
- (4) Inspiratory resistance  $\leq 175\text{ Pa}$  and expiratory resistance  $\leq 145\text{ Pa}$ .
- (5) The visual field under the mask is  $\geq 60^\circ$ .



3 ply disposable non-medical face mask with elastic




Main material: This product consists of a protective body, a nose clip and an ear strap. The protective body is made of nonwovens (three layers of the inner and outer layers and the middle filter layer)

Inner Layer -  
spunbond  
nonwoven fabric

Filter layer - melt-  
blown nonwoven

Outer layer -  
spunbond  
nonwoven fabric

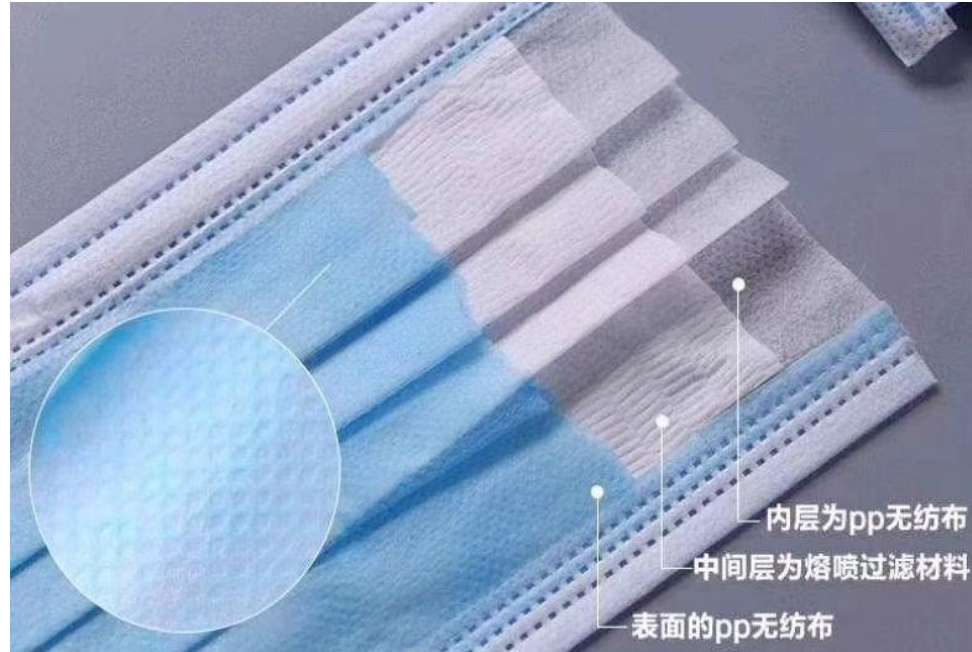

## Physical sample size





# Pacific Import Manufacturing, Inc.

**Calvin Wang**  
**President**

13000 Temple Avenue,  
City of Industry, CA 91746

**Phone: (818) 404-0806**  
**Email: [snwangpilot@earthlink.net](mailto:snwangpilot@earthlink.net)**